

# Quantum® 1107

POWER CHAIRS

Quantum® 1107 shown with  
Synergy® Manual Recline and  
PG Remote Plus controller



## Small frame. Big torque.

A highly stable and compact design gives the Quantum® 1107 outstanding rehab capability, while a powerful drivetrain delivers superb performance indoors and out. The Quantum 1107 is compatible with an impressive array of electronics and seating and positioning systems.



## Specifications

<b>DRIVE WHEELS</b>	10"
<b>REAR CASTERS:</b>	6" solid
<b>ANTI-TIPS:</b>	5"
<b>SUSPENSION</b>	Advanced Active-Trac® with rear caster suspension
<b>MAXIMUM SPEED</b> <sup>1</sup>	Up to 5 mph
<b>BRAKES</b>	Intelligent Braking (electronic, regenerative disc brakes)
<b>GROUND CLEARANCE</b>	2.5"
<b>TURNING RADIUS</b> <sup>3</sup>	21.875"
<b>OVERALL SIZE</b> <sup>3</sup>	Length: 36.375" Width: 23.375"
<b>SEATING</b>	Synergy® Seat Adjustable Width: 10"-20" Adjustable Depth: 10"-20" TRU-Balance® 2 Power Positioning System <sup>2</sup> Adjustable Width: 14"-20" Adjustable Depth: 14"-20" Specialty Seat
<b>DRIVE TRAIN</b>	2 motor, mid-wheel drive;
<b>BATTERIES</b>	Two 12 volt, deep cycle U-1
<b>RANGE PER CHARGE</b> <sup>1</sup>	Up to 15 miles
<b>BATTERY CHARGER</b>	5A, off-board
<b>AVAILABLE ELECTRONICS</b>	75A, Flight 70A, PG VSI 70A, Q-Logic
<b>WEIGHT CAPACITY</b>	300 lbs. maximum (user weight)
<b>BASE WEIGHT</b>	112 lbs.
<b>BATTERY WEIGHT</b>	24.5 lbs. (each)
<b>WARRANTY</b>	Lifetime limited warranty on frame 2-year warranty on electronics 18-month warranty on drive motor


## Shroud Colors



Candy  
Apple Red



Viper  
Blue

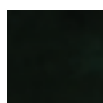

Onyx  
Black



Sunburst  
Orange

## Features

- GLIDE motors for added power and torque
- Innovative battery boxes and front module eliminate battery and motor harnesses
- Advanced Active-Trac and rear caster suspension for enhanced performance on rough and uneven terrain, creating a smooth and stable ride
- Standard 5 amp, Lester Power Cube off-board charger

## Options

- TRU-Balance® 2 Power Positioning System
- Q-Logic drive control system
- Contoured high-back seat
- Solid tires
- Cup holder
- Weather cover
- Swing-away joystick
- Elevating leg rests
- TRU-Comfort Seating
- Oxygen tank holder
- Cane/crutch holder
- Walker holder
- Rear basket
- Angle-adjustable footplates
- 70° heavy-duty drop-in footrests

## Code Coverage

Quantum 1107 2SP - SS | (Code: K0835)

(US) 866-800-2002  
(Canada) 888-570-1113  
www.quantumrehab.com



### FOOTNOTES:

- 1) Range and speed vary with user weight, terrain type, battery charge, battery condition and tire pressure.
- 2) 200 lbs. weight capacity with TRU-Balance 2 Power Positioning System
- 3) Due to manufacturing tolerances and continual product improvement, this specification can be subject to variance of + or - 3 %.